



Celesta



Shown with Clay top, Walnut cabinet and optional ABC System

The Celesta, the ultimate in luxury with exquisite design, superior comfort and fully electric functionality will transform the spa experience. Featuring Oakworks unique 4 electric motor system and now standard with a removable and replaceable padded top. The clean lines foster a sense of elegant simplicity, while the custom-designed cabinetry, complete with electrical outlets inside the cabinet, adds a functionality that makes incorporating a hot towel cabi a convenient reality.




Bluetooth Speaker System (Optional)

Delivers excellent quality stereo sound via any Bluetooth connected device.



Programmable Hand Control

This control has two user programmable memory positions and the ability to change two settings at once with a single button. The simple interface makes positioning clients fast and the new electronics are quiet and smooth in operation.

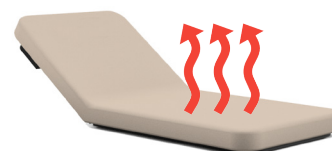

Integrated Heated Top (Optional)

Our new system is imbedded in the foam, creating a more even heat, and has three settings. The heating element automatically shuts off after 10 hours providing peace of mind for safety and energy conservation.

SPECIFICATIONS

- New! Removable and replaceable padding
- Programmable memory positions
- AeroCel™ padding for ultimate comfort
- Large height range - 25" - 35"

TABLE SPECIFICATIONS

Top	Electric Salon Top or Flat top
Lift System	4 Electric Motors
Controls	Foot Control + Hand Control with memory positions
Lifting Capacity	500 lbs. (249 kg)
Fabric	TerraTouch - Medical Grade, ISO 10993
Finish	Water Based Stain or Lacquer
Padding	4" (10 cm) or 6" (15 cm) AeroCel™ (Optional)
Height Range with 4" (10.2 cm) pad	25" - 35" (63 cm - 89 cm)
Width	29" (74 cm), 31" (79 cm), 33" (84 cm)
Length	73" (185 cm)
Options (must be ordered with table)	Integrated Heated Top, ABC System, Oakworks Bluetooth Speaker
Warranty	3 years parts, 2 years labor
Safety Listings	cETLus listed 3034177

Accessories	Adjustable Side Arm Rests, Boiance Shelf Support, Boiance Face Rest Cushion, Quick Lock Face Rest
-------------	---

STAINS



LACQUERS

TERRATOUGH COLORS

Meets the requirements of ISO 10993-1:2009 Biological Evaluation of Medical Tables & passes testing for Cytotoxicity, Skin Irritation & Sensitization as well as CA TB117 Fire Retardance testing. Please refer to the product cleaning instructions for a list of approved cleaners and disinfectants.



ABC SYSTEM (OPTIONAL)

Oakworks Patented ABC system allows you to create the perfect level of support for each client, taking pressure off their breasts and tension out of their back. The system also raises the pad above level, making it a great head support for face up positions, and can be used for stretches also.

(Comes with Boiance face rest pad and QuickLock face rest)



Incredible breast comfort
Lower to lessen pressure



Electric adjustments
Large adjustment range



Head support
Built in head support



Pectoral muscle stretch
Advanced positioning

ADJUSTABLE SIDE ARMREST (ACCESSORY)

These versatile and comfortable removable arm rests feature a wide range of positions. They always stay level, regardless of the backrest angle. They offer comfortable positioning and can be moved completely out of the way to allow guests to get on and off the table easily.



BOIANCE SHELF SUPPORT (ACCESSORY)

The shelf and bracket used together allow you to create the perfect place for clients arms when working on their back. The shelf can also be used alone as a head support for facials and other aesthetic treatments. The shelf is designed to allow therapists to relax their shoulders while working on clients, minimizing stress and injuries.



Contact Oakworks
717.235.6807
www.shop.spatable.com
customerservice@oakworks.com