


## PF SERIES EXAM TABLES

The PF Series offers a great combination of comfort and beauty for non-invasive procedures in high-end treatment spaces. The tilt function allows you to create a Zero Gravity experience that is enhanced by the plush foam system. With it's beautiful wood construction and different finishes you can create a Spa like experience in a medical environment.



PF400 shown with sapphire fabric and optional grey lacquer

### One Piece Replaceable Pad

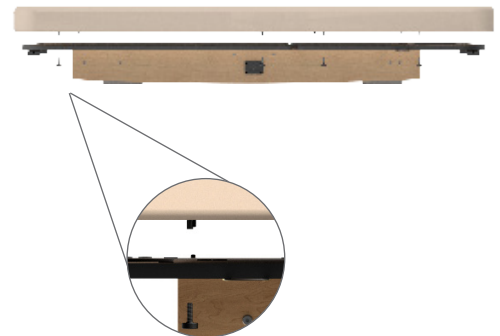

Our design features finger knobs that can be removed in minutes, allowing you to completely remove the entire top and install a new one without any special tools so your table will be in service for much longer.



### Fully Electric Contour Top



## MODEL PF400

Motion 1	Electric Height: 21"-37" (53-97cm)
Motion 2	Electric Tilt: +/-15°
Motion 3	Electric Backrest: 0-55°
Motion 4	Electric Leg rest

### Lift-Assist Backrest Top



## MODEL PF250

Motion 1	Electric Height: 22"-37" (56-95cm)
Motion 2	Electric Tilt: +/-15°
Motion 3	Lift Assist Backrest 0-80°

# SPECIFICATIONS


- Creates a Spa experience in your practice
- Tiltable top with luxurious padding for greater comfort
- Large height range - 21" - 37" (3" pad)
- Beautiful wood construction with choice of water resistant finishes

## TABLE SPECIFICATIONS

Lift System	Electric, Twin Tower
Controls	Foot Control, Hand Control (Optional)
Weight Capacity	550 lbs. (249.5 kg)
Fabric	TerraTouch
Finish	6 Water Resistant Finishes on Wood
Padding	3" (7.6 cm) or 4" (10.2 cm) Comfort Foam (Optional)
Motion 1 - Height Range 3" (7.6 cm) Pad	PF400 - 21" - 37" (53 cm - 94 cm) PF250 - 22" - 37" (56 cm - 95 cm)
Motion 2 - Tilt	+/-15°
Motion 3 - Backrest	PF400 - Electric Backrest: 0-55° PF250 - Lift Assist Backrest: 0-80°
Motion 4 - Leg Rest	PF400 - Electric Leg Rest
Width	29" (73.6 cm), 31" (78.7 cm), 33" (83.8 cm)
Length	73" (185.4 cm)
Options (must be ordered with table)	Integrated Heated Top, Casters, Hand Control
Warranty	3 years parts, 2 years labor
Safety Listings	FDA listed, cETLus listed 3034177, CB Certificate


Accessories	Adjustable Side Arm Rests, Quicklock Head Rest & Aero Cel Pad, Universal Table Extender
-------------	---

## WOOD FINISHES

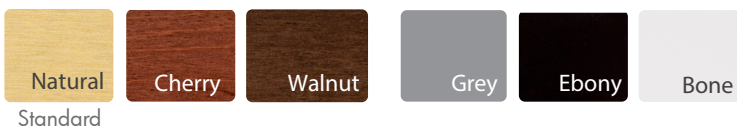


## TERRATOUGH FABRIC COLORS

Meets the requirements of ISO 10993-1:2009 Biological Evaluation of Medical Tables & passes testing. Please refer to the product cleaning instructions for a list of approved cleaners and disinfectants.

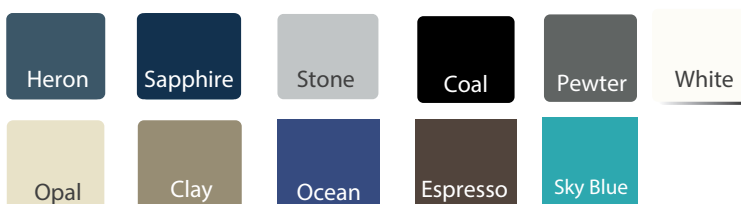


## TRENDELENBURG TILT (STANDARD)

The unique twin tower lift system creates Trendelenburg tilt capability insuring the ultimate in comfort. The tilting function is helpful for many treatments and offers greater comfort for people who experience sinus pressure when lying flat.

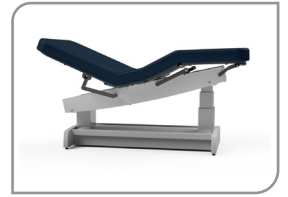

## INTEGRATED HEATED TOP (OPTIONAL)

Our new system is imbedded in the foam, creating a more even heat, and has three settings. Designed originally to be used in car seats and hospital beds, the heating element automatically shuts off after 10 hours providing peace of mind for safety and energy conservation.

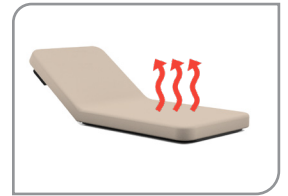

## 3" LOCKING CASTERS (OPTIONAL)

Locking casters make it easy to move the table in the room.



## QUICKLOCK™ HEAD REST & AERO·CEL PAD (ACCESSORY)

The platform allows precise cervical flexion and height positioning over a wide range for both prone & supine treatments.




## UNIVERSAL TABLE EXTENDER (ACCESSORY)

Perfect as a head rest or a foot rest, it has two rods that fit into the bracket on the table. Extends table to 85".



## ADJUSTABLE SIDE ARMREST (ACCESSORY)

These versatile and comfortable removable arm rests feature a wide range of positions. They always stay level, regardless of the backrest angle. They offer comfortable positioning and can be moved completely out of the way to allow patients to get on and off the table easily. (PF400 only)

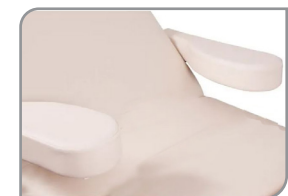

717.759.3146 [www.oakworksmed.com](http://www.oakworksmed.com) email: [medsales@oakworks.com](mailto:medsales@oakworks.com)